

**The Leaf Fusion Plasma Device** features 6 main skin care functions and technologies in one single device. Use it in the comfort of your own home.

This unique, patented, and state-of-the-art skin care solution is based on the following three operating modes (each individual mode with three intensity levels):

FUSION PLASMA MODE	SMART BOOSTING MODE	MICRO-CURRENT MODE
<ul style="list-style-type: none"><li>• Securing stratum corneum passage</li><li>• Wrinkle improvement</li><li>• Rough skin texture improvement</li></ul>	<ul style="list-style-type: none"><li>• Accelerate skin care product application</li><li>• World's first, 3rd generation electroporation technology</li></ul>	<ul style="list-style-type: none"><li>• Lifting</li><li>• Promoting cell regeneration</li></ul>


In addition, the following three technologies are also incorporated into the device:

VIBRATION	PROBE WARMING	PHOTO FACIAL
<ul style="list-style-type: none"><li>• Skin muscle relaxation</li><li>• Blood flow promotion</li></ul>	<ul style="list-style-type: none"><li>• Reserve collagen</li><li>• Promote blood flow and metabolism</li></ul>	<ul style="list-style-type: none"><li>• LED simultaneous operation</li><li>• Promoting blood flow and anti-aging</li></ul>



Plasma Arc + Electroporation  
Ozone Free

Smart Update Function  
Convenient voice guidance function  
Dual charging by Cradle/ Main Body

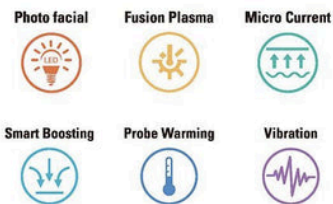

## Contents

Leaf Fusion Plasma device

Charging cradle

Micro USB charging cable

Manual and warranty card

1.65 (W) x 7.48 (H) x 1.57 (D) inches = 42(W) x 190(H) x 40(D) mm

## Smart Update/Calibrate Function via Mobile Phone

Select the best course according to the type of skin care product by mobile phone screen communication.

By updating the information with your specific skin care product and skin type, you can customize it and increase the absorption rate



## Skin Sensing Function

The LED flashes and vibrates only when contact is made with the skin. If the skin is too dry, the device may not work properly. Please apply enough skin care product. When proper skin condition is sensed, the LED will flicker & vibration starts




## Unique patented skin care solution in three steps:

Exfoliate and promote absorption of skin care products through fusion plasma and electroporation. Face lifting through micro-current massage.

### 1. FORMATION OF THE STRATUM CORNEUM THROUGH FUSION PLASMA


Fusion plasma energy creates a microscopic pathway in the stratum corneum layer of the skin. Since the passageways are formed at a microscopic depth, a primary passage is secured to absorb the active ingredients of the cosmetic without damaging the skin.

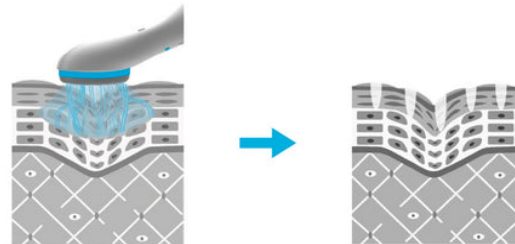

It also stimulates the external stimulation of the skin by heat of the fusion plasma, induces skin regeneration and promotes skin lifting.


### 2. SECURING PASSAGE BY 3RD GENERATION ELECTROPORATION AND INPUTTING ACTIVE INGREDIENT

Using 3rd generation electroporation, it maintains the passage formed in the first step and secures additional passage to promote absorption of active ingredient in cosmetics. The third-generation electroporation primarily increases the permeability of the field strength and cell membrane, and secondarily moves the active ingredients of the cosmetic product down the stratum corneum to promote absorption.

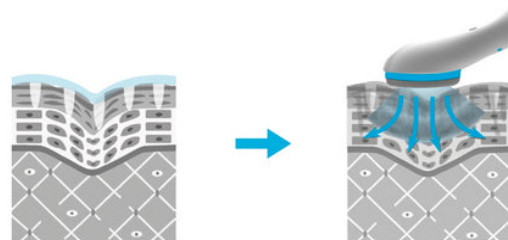

### 3. LIFTING USING MICRO-CURRENT

It is a step that stimulates lifting by helping to regenerate ingredients related to skin elasticity by applying lower electric stimulus than the previous methods. Collagen and elastin-related ingredients, such as elasticity, age begins to slow down the recovery. At this time, the energy of the micro-current delivered to the skin through electrical stimulation stimulates the ingredients to help regenerate and promote skin lifting.

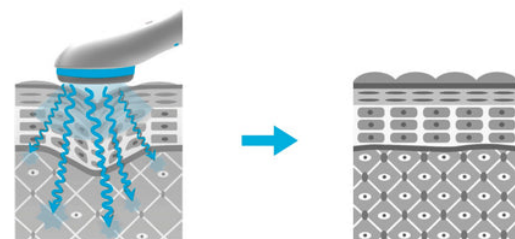

## **HOW TO USE**

### **1. FUSION PLASMA STEP**

1. Fusion plasma step to secure stratum corneum passage:
2. Facial wash with cleansing cream or foam.
3. With nothing applied wait until the water dries slightly.
4. After turning on the power[ ], it begins from 1st intensity. To select Fusion Plasma mode in other modes, press shortly to change to Fusion Plasma mode.
5. As the probe hardly touches the surface of the skin, lightly massage the entire face as if brushing.

Please note: If there is water on your face, plasma does not form properly. Please dry thoroughly before proceeding with Fusion plasma step.

### **2. SMART BOOSTING STEP**

1. Promoting cosmetic application by skin boosting by Electroporation:
2. Apply proper amount of cosmetics evenly over the face.
3. Press smart boosting mode button. It begins from 1st intensity.
4. Put the probe in close contact with the face and evenly boost whole face as if massaging for about 5 mins. The upper part of the head makes it easier to care for narrow areas.
5. You can adjust the use time according to the degree of absorption of cosmetics.
6. Depending on the amount of oil and moisture in your skin, it may stop working. Apply sufficient amount of cosmetics and make the optimal condition for absorbing cosmetics.
7. When using on the neck, please be careful not to directly affect the uvula.

### **3. MICRO-CURRENT LIFTING STEP**

Lifting by micro-current:

1. Facial wash with cleansing cream or foam.
2. Apply proper amount of cosmetics for lifting and anti-aging evenly over the face.
3. Press shortly the lifting mode button. It begins from 1st intensity.
4. Apply the head of the device to the surface of the face for 5 minutes and slowly massage it while lifting it along the direction of the arrow as in the picture.
5. You can adjust the use time according to the degree of absorption of cosmetics.